




Accel® TB

Infection Prevention with Confidence

In the war against microbes found on surfaces in critical environments and controlled areas, Accel® Disinfectants are the first line of defense for infection prevention and control.

EPA registered Accel TB Disinfectants, offered exclusively by Contec, are available in ready-to-use solutions and presaturated wipes for easy disinfection of environmental surfaces and patient care devices. Accel TB Products are powerful intermediate disinfectants that utilize Accelerated Hydrogen Peroxide® (AHP®) technology and are effective against a broad-spectrum of pathogens in only 1 minute!

Accel TB is ideal for use on environmental surfaces in pharmaceutical cleanrooms, compounding pharmacies, labs, and many other GMP facilities. The ready-to-use formula is also convenient for forensic/crime scene clean up, and other controlled environments.

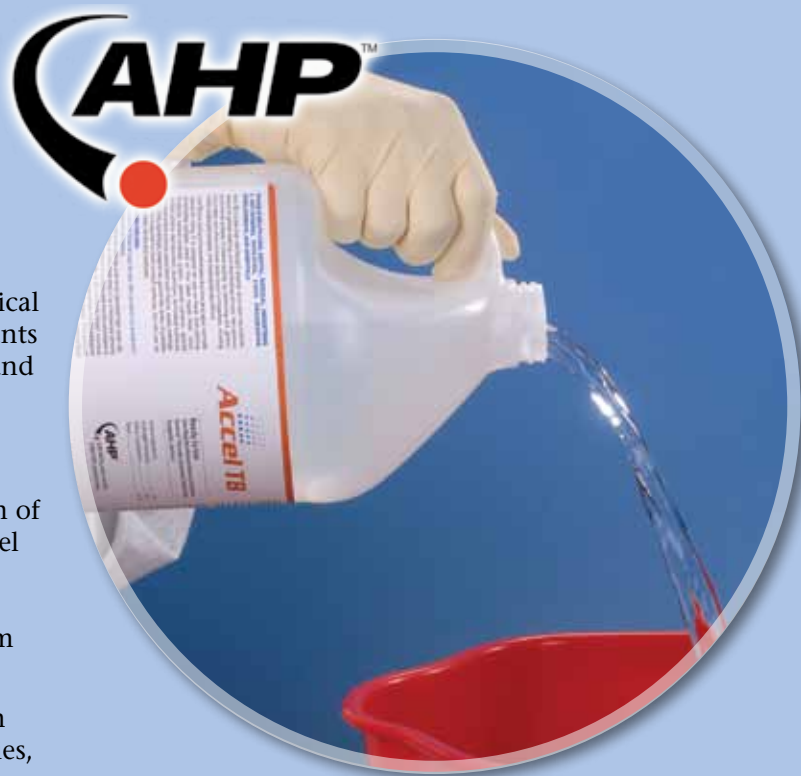

Accel® TB from Contec® is made with AHP, Accelerated Hydrogen Peroxide, a patented cleaning and disinfecting solution.



EPA Registered Claims

- 30 Second Non Food Contact Broad-Spectrum Sanitizer giving a 3 Log (99.9%) reduction in Vegetative Bacteria*
- 1 minute Bactericidal Claim giving a 6 Log (99.9999%) reduction in Vegetative Bacteria*
- 1 minute General Virucide Claim giving a 3 Log (99.9%) reduction*
- 5 minute Mycobactericidal Claim giving a 4 Log (99.99%) reduction*
- 10 minute Fungicidal Claim giving a 5 Log (99.999%) reduction*

**In the presence of 5% serum load at 20°C on hard, non-porous environmental surface. Please see technical data sheet or website for a complete list of organisms and viruses.*





For more documents about Accel TB, please visit the Resources section on our website to download:

- AHP Brochure
- Technical Data Sheet
- Material Safety Data Sheet
- Compatibility Chart
- Press Release and Articles
- Product Labels



P.O. Box 530 • Spartanburg, SC 29304 USA
+1-864-503-8333 • fax: +1-864-503-8444
Toll Free: 1-800-289-5762
www.contecinc.com • e-mail: disinfect@contecinc.com

Features & Benefits

FEATURE	BENEFIT
 <p>FASTER <i>Realistic Contact Times</i></p>	<ul style="list-style-type: none"> • Short contact times on a broader spectrum of organisms ensures user compliance.
 <p>CLEANER <i>Better Cleaning Efficiency</i></p>	<ul style="list-style-type: none"> • Better cleaning efficiency results in added confidence that disinfection can occur. • No active residues left on the surface as the active ingredient breaks down into water and oxygen
 <p>SAFER <i>Safer for Users and Occupants</i></p>	<ul style="list-style-type: none"> • All chemicals are on the FDA Generally Regarded as Safe List (GRAS). • Chemicals that are safer for users provide for a healthier indoor environment.
 <p>RESPONSIBLE <i>Reduced Environmental Load</i></p>	<ul style="list-style-type: none"> • Reducing the load on the environment whenever possible is a responsible thing to do. • Sustainable solutions contain safer ingredient choices.

Accel® TB Disinfectant Solution

Ready-to-use liquid with 0.5% Hydrogen Peroxide, EPA # 74559-1

Accel TB Disinfectant Solution is an EPA-registered intermediate level disinfectant that is safe and fast acting with a broad spectrum kill. Accel TB utilizes AHP®, a synergistic blend of commonly used, safe ingredients that, when combined with low levels of hydrogen peroxide, dramatically increase its germicidal potency and cleaning performance. Compatible with an extensive list of common materials, including 316 and 304 stainless steel, natural and silicone rubber, and most plastics, Accel TB is ready-to-use straight from the bottle, requires no dilution, and has a 2-year shelf life; opened or unopened.



Part Number	Description	Packaging
ACCDISR-TB32	Accel TB Ready to Use Liquid, 32 oz. (0.9L) bottle	12 bottles/case
ACCTRIGR	Non-aerosol trigger sprayer for 32 oz. bottle	12 sprayers/case
ACCDISR-TBG	Accel TB Ready to Use Liquid, 1 gallon (3.8L) container	4 containers/case
ACCDISR-TB5G	Accel TB Ready to Use Liquid, 5 gallons (18.9L) container	1 container/case

Accel® TB Disinfectant Wipes

Presaturated wipes with 0.5% Hydrogen Peroxide, EPA # 74559-3

Exclusively available through Contec, Accel TB Disinfectant wipes combine the same efficacy as Accel TB solution with the convenience of presaturated wipes. Ideal for use in areas with limited storage or counter space, Accel TB wipes measure 6" x 7" (15.2 x 17.8cm), are packaged 160 wipes per canister, and have a 3-year shelf life. Presaturated wipes are made with clean processed meltblown polypropylene and have the added benefit of standardizing work habits; the same amount of disinfectant solution is released with each wipe.



Part Number	Description	Packaging
ACCWIP1-TBUS	Accel TB Disinfectant Wipes, 6" x 7" (15.2 x 17.8cm)	160 wipes/canister 12 canisters/case